

# Jazzy<sup>®</sup>

SELECT 6 ULTRA

## Jazzy<sup>®</sup> Select<sup>™</sup> 6 Ultra

The Jazzy<sup>®</sup> Select<sup>™</sup> 6 Ultra provides the advantages of in-line motor technology, maximum front stability and the environmental access benefits of a standard power elevating seat. These features, along with Active-Trac<sup>®</sup> Suspension and 6" front and rear casters make the Jazzy Select 6 Ultra high performing and easy to use.



JAZZY<sup>®</sup> SELECT<sup>™</sup> 6 ULTRA shown in Jazzy Red with contemporary high-back seat, PG 6-key VR2 controller and foot platform



The Jazzy® Select™ 6 Ultra series models include:

- Jazzy Select 6 Ultra-C
- Jazzy Select 6 Ultra-SS

#### FEATURES

- Standard power elevating seat
- 6" casters on front & rear for maximum stability
- Dual in-line motors for enhanced efficiency, torque, range and performance
- Built with ease of service in mind
- Battery access through side panel without the removal of the seat
- Tool-less front panel access to electronics to minimize service time
- New 40 amp, PG 6-key VR2 controller

#### OPTIONS

- High-back seat with solid seat pan
- Cup holder
- Weather cover
- Rear basket
- Oxygen tank holder
- Cane/crutch holder
- Walker holder

The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3 %.

Jazzy® Select™ 6 Ultra	
Weight capacity:	300 lbs.
Maximum speed:	Up to 4 mph
Ground clearance:	2"
Turning radius:	20.5"
Overall length:	35.625" without front riggings
Overall width:	23.5"
Wheels:	Drive: 10"
Front casters:	6" solid
Rear casters:	6" solid
Suspension:	Active-Trac®
Drive train:	Two-motor, in-line, mid-wheel drive
Braking system:	Intelligent braking (electronic, regenerative disc brakes)
Std. Electronics:	40 amp, PG VR2 controller
Specialty controls:	N/A
Battery charger:	O -board, 5 amp
Per-charge range:	Up to 15 miles
Battery requirements:	(2) 12 volt, deep cycle
Battery size:	U-1
Transport weight	
Battery weight:	24.5 lbs. each
Base weight:	114 lbs.
Max. Seat dimensions:	20" x 20"
Warranty:	5-year limited warranty on frame; 1-year limited warranty on electronics; 1-year limited warranty on drive motors

Colors:




#### Pride Silver Star Lifts & Ramps

- Full Platform
- Backpacker® Plus
- Powerchair Lift with straps



Backpacker® Plus

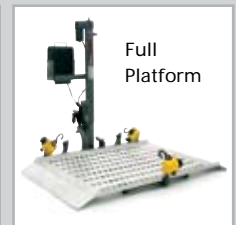

Full Platform

These lifts are compatible with the Jazzy® Select™ 6 Ultra.