

GO-GO[®]

LX WITH CTS SUSPENSION

3 OR 4-WHEEL SCOOTER




Comfort-Trac Suspension

Features, Performance and Travel Ready

The Go-Go[®] LX with CTS Suspension (Comfort-Trac Suspension) brings a new level of performance, features and value to travel mobility. Featuring CTS independent front and rear suspension, sleek style, feather-touch disassembly and standard lighting, the Go-Go LX with CTS Suspension is the most advanced travel product available.



Pride[®]
Mobility Products Corp.
— LIVE YOUR BEST —

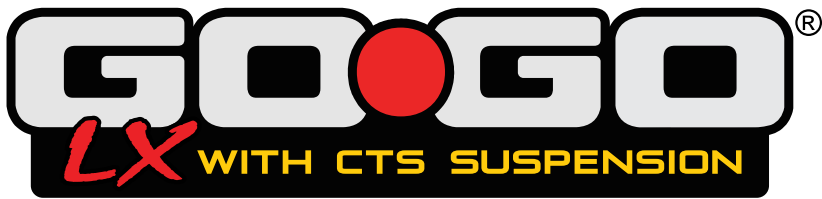

SPECIFICATIONS

Model number	S50LX (3W), S54LX (4W)
Max. weight capacity	300 lbs.
Maximum speed ¹	Up to 4.4 mph
Ground clearance ²	1.5" at motor
Turning radius ²	33" (3W), 44.75" (4W)
Overall length ²	39.75" (3W), 41" (4W)
Overall width ²	20.5"

TIRES:

Front	2" x 8" solid (3W), 2" x 7" solid (4W)
Rear	2.5" x 8" solid
Anti-tips	Rear

Standard seating	Foldable, compact, vinyl-covered molded plastic
------------------	---

Drivetrain	24-volt DC motor, sealed transaxle
------------	------------------------------------

Dual braking system	Regenerative and electromechanical
---------------------	------------------------------------

Battery charger	Off-board, 2 amp
-----------------	------------------

Per-charge range ^{1,4}	Up to 9.44 miles (3W), Up to 10.5 miles (4W)
---------------------------------	--

Battery requirements ³	(2) 12V deep-cycle
-----------------------------------	--------------------

Compatible batteries	(2) 12AH (standard) or (2) 18AH
----------------------	---------------------------------

WEIGHT:

Batteries (2 required) ⁶	12AH (9 lbs. each), 18AH (14 lbs. each)
Weight w/o battery pack	85.4 lbs. (3W), 87.4 lbs. (4W)
Heaviest piece	34 lbs. (front section, 3W) 36.6 lbs. (front section, 4W)

Interchangeable shroud panels



Red

Blue

FEATURES

- CTS Suspension (Comfort-Trac front and rear independent suspension)
- Feather-touch disassembly permits simple frame separation with only one hand
- Pride's exclusive black, non-scuffing tires
- Auto-connecting front-to-rear cable
- Includes 2 sets of easily changeable colored panels in Red and Blue
- Front frame-mounted seat post offers maximum stability
- All-new frame design easily disassembles into 5 super lightweight pieces for convenient transport and storage
- Modular design for easy serviceability
- Microprocessor-based controller offers optimal power management and added safety features
- Standard front tiller-mounted basket included
- Convenient off-board dual voltage charger can charge battery pack on-board or off-board
- Charger port in tiller

Pride® Lifts & Ramps

These Pride® Lifts are compatible with the Go-Go® LX with CTS Suspension.

- Outlander Exterior Lift
- Power Tote Mini



Power Tote Mini



Outlander

¹Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tire condition.

This specification can be subject to a variance of (+ or -) 10% - 5%.

²Due to manufacturing tolerances and continued product improvement, this specification is subject to a variance of (+ or -) 3%.

³AGM or gel cell type required.

⁴Tested in accordance with ANSI/RENSA, WC Vol. 2, Section 4, and ISO 7176-4 standards. Results derived from theoretical calculations based on battery specifications and drive system performance. Test conducted at maximum weight capacity.

⁵Includes standard seat.

⁶Battery weight may vary based on manufacturer.

The information contained herein is correct at the time of publication. We reserve the right to alter specifications without prior notice.