



Shown in Maple

## Standard Fabrics:



Acorn



Maple

This faux leather is a rich and modern alternative to natural leather that is extremely comfortable while offering the rugged durability of leather. Made from 100% polyester, the fabric is easy to care for ensuring functionality as well as elegance. This long-lasting fabric will complement a variety of decors while standing the test of time under heavy use.

Crossover fabrics are unavailable.



**Rhea  
PR442**



*Comfort Zone™*  
by GOLDEN

Take a vacation from your busy day in The Rhea. It offers an overstuffed seam back design, plush chaise pad and an open arm design that holds numerous standard features that include two cup holders, two storage compartments, and a USB port. Positioning options are endless with our 4-Zone Comfort Zone™ Positioning System, which includes our patented MaxiComfort® Positioning, adjustable headrest and lumbar all operated by our programmable AutoDrive™ hand control. For extra wellness and comfort, our HeatWave™ Technology comes standard providing therapeutic warmth starting at your shoulders and continuing down to your feet. When you open the right arm, you will find a cup holder and a storage compartment that also holds the AutoDrive hand control. In the left arm, you will find an additional cup holder, a USB charging port, and a second storage compartment that holds the HeatWave hand control. As a bonus, you can use either cup holder without losing the comfort of the armrests thanks to the split design!

## Standard Features:

- Patented MaxiComfort® Positioning
- Adjustable Headrest
- Adjustable Lumbar
- Programmable AutoDrive™ Hand Control
- HeatWave™ Technology
- USB Charging Port
- Built-in Cup Holders
- SmartTek® Diagnostic System
- Extra Pocket



# DeLUNA Series PR442

Model Number	PR442-MED	PR442-LAR
Size	Medium	Large
Class II Medical Device	YES	YES
Weight Capacity	375 LBS	375 LBS
Positioning	MaxiComfort	MaxiComfort
Comfort Zone(s)	4	4
Twilight Technology	NO	NO
Chaise Seat	YES	YES
Pocketed Coil Seat	YES	YES
Back Type	Seam	Seam
Removable Back Type	Bracket	Bracket
Head Rest Cover	NO	NO
Arm Rest Covers	NO	NO
Overall Width	37"	39"
Overall Height	45"	50"
Width Between Arms	22"	24"
Floor to Top of Seat	20"	21"
Seat Depth	21"	22"
Seat to Top of Back	29"	33"
Seat to Top of Arm	7"	8"
Distance Required from Wall Reclined	25"	30"
Extension Reclined from Back of Seat to End of Foot Rest	43"	44"
Life Time Warranty, Chair Frame, Lift Frame, Recline Mechanism	YES	YES
Warranty on Electrical Parts	3 YR	3 YR
Pro Rated Years 4 Thru 7	YES	YES
Warranty on Mechanical Labor	3 YR	3 YR
Warranty on Electrical Labor	1 YR	1 YR
Hand Control Routing	Right Arm Compartment	Right Arm Compartment
Extra Pocket Option	Standard	Standard
HeatWave Technology	Standard	Standard
Heat & Vibrating Massage Option	NO	NO
Articulating Headrest Option	Standard	Standard
Adjustable Lumbar Option	Standard	Standard
AutoFlex Option	NO	NO
Foot Rest Extension Option	NO	NO
Straight Lift Option	YES	YES
USB Charging Port	Left Arm Compartment	Left Arm Compartment
Battery Backup	Lithium-Ion	Lithium-Ion
Place of Origin	China	China

SmartTek Diagnostic System STANDARD on all Golden models.

Designed and engineered in Old Forge, PA. Assembled in China using Golden components.



**HEAT WAVE**  
TECHNOLOGY

### STANDARD!


Far infrared heat penetrates skin and muscles to provide a soothing heat therapy experience. Adjustable Hi-Lo settings and timers let you control the comfort!



The Rhea comes standard with the following comfort and convenience features:



Right Arm: Cup Holder & Storage Compartment

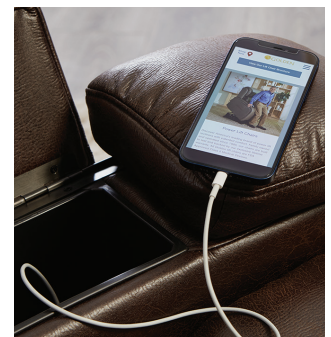

Left Arm: USB Charging Port, Cup Holder & Storage Compartment



Adjustable Headrest & Lumbar



HeatWave™ Technology




## Comfort Zone™ by GOLDEN

Comfort Zones 1 - 4 provide an ergonomic comfort experience with an emphasis on head and back support. Each Comfort Zone can be adjusted independently by using the exclusive AutoDrive hand control.



LOWER BODY UPPER BODY HEADREST LUMBAR TWILIGHT TILT

**WE PATENTED COMFORT™**

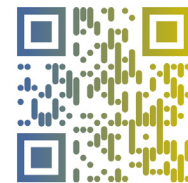

Scan with your smartphone camera to learn more about Comfort Zone by Golden.



Phone: 800-624-6374  
Fax: 800-628-5165  
www.goldentech.com

Literature is current at time of printing. Golden Technologies reserves the right to make changes to the product or literature at any time.

Sizes may vary depending on fabric, filling material, upholstery and carpet thickness. All measurements taken on concrete floor using levels and metal rulers.



Complies with  
UL 73 and  
CSA C22.2 No. 68



www.facebook.com/goldentechnologies

SHT-PR442-081623